
Universal Silberlotpaste
Universal Silver Solder Paste
Brasure d'argent universelle en pâte
Soldadura universal de plata en pasta
Saldatura universale in pasta

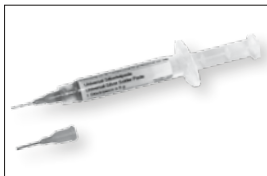Gebrauchsanweisung	S. 3
Instructions for use	P. 7
Mode d'emploi	P. 11
Modo de empleo	P. 15
Modalità d'uso	P. 19

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Dear Customer

Thank you for having chosen a quality product from Dentaureum.

In order to use this product at its best for you and your patients, it is important to study and follow these directions for use carefully.

The written instructions cannot cover all eventualities during operation. For this reason our Hotline is available to answer any other questions and ideas that may arise.

Due to a constant development we recommend, even when you use the same products frequently, to study the enclosed latest issue of the instructions for use or refer to our website at www.dentaureum.de.

Manufacturer

Dentaureum J. P. Winkelstroeter KG
Turnstraße 31 · 75228 Ispringen · Germany

With reference to our quality standards

Dentaureum ensures a faultless quality of the products manufactured by us. These recommendations are based upon our own experiences. The user himself is responsible for the processing of the products. Responsibility for failures cannot be taken, as we have no influence on the processing on site.

Instructions for use

Universal Silver Solder Paste

This silver solder paste has been specially developed for use in orthodontics.

Its distinct feature is its easy and precise application and bonus of requiring minimal quantities. The flux which is integrated within the paste is exactly suited to the solder, so that no extra flux is required.

When using the universal silver solder paste, excellent results can be achieved by both the torch soldering and electro soldering techniques.


The comparatively low working temperature (630 °C / 1166 °F) means it is easy to handle.

1. Application

All types of commercially standard stainless steel can be soldered with this paste.

2. Working instructions

- a) Remove the lid and replace with the nozzle.
- b) Apply the universal silver solder paste onto the necessary areas on the piece of work (solder gaps) by carefully pressing the syringe plunger into the cylinder.
- c) When the paste has been applied to the areas to be soldered, heat the piece of work using a spot welding unit (e.g. REF 079-000-00) or a soldering torch (working temp. 630 °C/1166 °F) until the solder flows and melts.
- d) Afterwards, if necessary, polish in the usual manner.

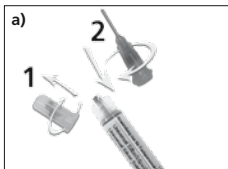3. Tips

- there is enough flux contained within the universal silver solder paste. It is not necessary to add more
- the area to be soldered must, as always, be free from grease, oxides and dirt
- the universal silver solder paste must not be applied to the piece of work during heating as it is possible that the nozzle may become blocked.

4. chemical components

in weight %: Ag 54 / Cu 22 / Zn 17 / Sn 7

5. Technical data

Melting interval	620 °C – 650 °C / 1148 °F – 1202 °F
Flux – and binding agent	54 %
Density of the solder alloy	9.5 g/cm ³
Color	greyish brown

6. Storage

36 months

7. Availability

	quantity	REF
Universal silver solder paste in syringe with flux, cadmium-free	8 g	380-804-50

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